


The Ferno XTRIC8 Rescue Sled is a versatile and cost-effective solution for universal rescues including lowland, mountain, confined space and high angle extrications.

The XTRIC8 is a compact and durable system. Smart features decrease the time and effort necessary to package patients while providing traditional functionality and tactical evacuation methods in an easy to use device.

## **Key Features**

- XTRIC8 is supplied as a complete, ready to use unit with all necessary hardware
- The XTRIC8 is pre-rigged for immediate deployment out of the bag
- Includes an internal patient harness for maximum patient security during vertical lifts
- Colour coded, 4-point lifting bridle is pre-set to aid faster implementation
- Reflective external webbing for enhanced visibility in low light conditions
- Straps and buckles are glove-friendly and designed for efficient attachment and release
- Robust carry bag with reflective webbing, four large pockets and adjustable backpack harness
- XTRIC8 is intuitive to use and minimises user training

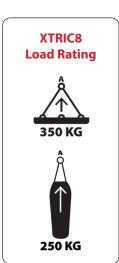





XTRIC8 RESCUE SLED	
Product Code	FWE XTRIC8
Inclusions	XTRIC8 Rescue Sled pre-rigged with internal patient harness assembly, external restraints, external straps and handles
	1 x Maillon Rapide for vertical lifting
	1 x Colour coded 4-Point Lifting Bridle
	1 x Storage Bag with backpack straps
Length	2540 mm
Width - Widest Point	914 mm
Weight - XTRIC8 only (rigged)	7 kg
Weight - XTRIC8 (rigged) and Bag	8.5 kg
Load Limit - XTRIC8	350 kg - HORIZONTAL / 250 kg - VERTICAL
XTRIC8 LIFTING BRIDLE	
Weight - Lifting Bridle	0.9 kg
Load Limit - Lifting Bridle	500 kg
Bridle Karabiners	VYO N-2853G - Aluminium Alloy (SAE 7075) Triple Action Karabiner
Bridle Karabiner Dimensions	124 mm x 80 mm / Gate Opening: 22 mm
Bridle Karabiner Ratings	Minimum Tensile Strength: 23 kN / Minimum Gate Strength: 16 kN
Bridle Karabiner Standards	ANSI Z359.12-2019; CSA Z259.12-16; EN 362:2004

• Ferno and manufacturers of components reserve the right to change specifications without notice.

• General specifications are rounded to the nearest whole number. Metric conversions are calculated before rounding the English measurements. For more information contact Ferno Customer Service or your local distributor.