




CapLiNXX

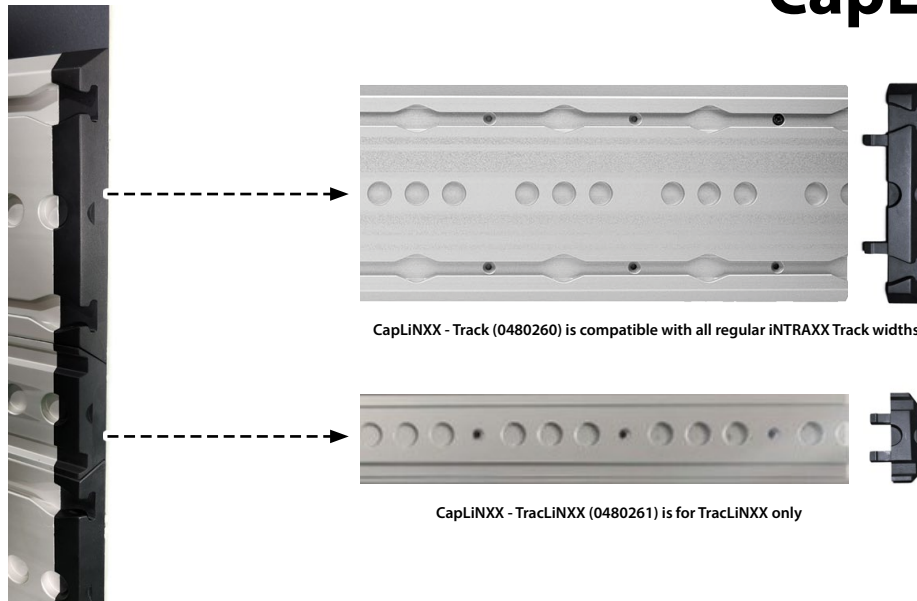

CapLiNXX, the newest addition to iNTRAXX's End Caps range, provides a functional, safe, and visually streamlined finish to the track. CapLiNXX enhances work efficiency and convenience when moving equipment from iNTRAXX Tracks.

The CapLiNXX profile allows mounts to be installed or removed from the end of the track groove without repositioning them to drop-in zones, all while maintaining a slim track profile. There are two CapLiNXX models —one designed for standard iNTRAXX tracks, and another for the TracLiNXX track profile.

Features

- Allows mounts to be installed or removed from track ends
- Lower profile is ideal for tight spaces
- Designed for reuse to reduce landfill waste, CapLiNXX secures to tracks without the need for additional fastening hardware
- CapLiNXX's low profile enables all tracks, regardless of length, to be 'capped' off without compromising their function — even in cases where one row of tracks is shorter than surrounding tracks
- Easy to clean and decontaminate
- Compatible with all iNTRAXX track types

Product Part #	0480260	0480261
Product Description:	CapLiNXX - Track	CapLiNXX - TracLiNXX
Overall Size: mm / inches	(L) 31 x (W) 189 x (D) 12.7 mm (L) 1.2 x (W) 7.4 x (D) 0.5 in.	(L) 31 x (W) 87 x (D) 12.7 mm (L) 1.2 x (W) 3.4 x (D) 0.5 in.
Weight	34 grams 0.1 lbs	16 grams 0.03 lbs
Load Capacity	N/A	N/A
Rating	10g, 22g, 26g	10g, 22g, 26g

